

MUSIC MEMORY

PRODUCT AVAILABILITY: Synergy Devices *i2400, i2000, and i1600 devices only

The Music Memory is a dedicated memory designed to enhance sound quality for patients who spend time listening to music. It complements the Auto Music feature to form Music Enhancement, a suite of features designed to improve music sound quality in a variety of environments.

Fitting and evaluation of the Music Memory are subjective. Make adjustments based on patient feedback regarding music sound quality. It is recommended to play music at the patient's preferred listening level while adjusting the Music Memory for best fitting results. Select the **Media Player** on the top tool bar, then **Music**, then select desired **Audio File**, or play music from patient's personal music device.

Devices feature a dedicated music memory, including:

-Compression architecture designed and optimized to enhance music sound quality


-Fitting rationale for music, focusing on providing audibility for soft music and transparent sound quality for loud music

-Inspire software adjustments tailored for music


1. Launch **Inspire X**, then connect and detect hearing devices. Refer to the *Initial Fitting Protocol QuickTIPs* for more information.
2. Select **QuickFit** on the left navigation bar.

MUSIC MEMORY

1. Select **memory 2, 3, and/or 4**, select the **memory drop-down arrow**, then select **Music** to program the Music Memory in one or more of the memory slots.
2. Select the **Music EQ Controls** on the QuickFit screen to adjust the Music Memory to the patient's preference. Adjustments include:
 - » **All Hz**
 - » **Bass**
 - » **Mid**
 - » **Treble**
3. Select **Fine Tuning** on the left navigation bar or from the QuickFit screen to access more channel-specific and input level-specific adjustments, if desired.
4. Reference Loud and Soft Music Equalizers to view Music Memory settings.
5. Select **+ Bass, Mid, and/or Treble** to expand channel ranges.
6. Select the **<< expansion arrow** to collapse the left navigation bar and enable full-screen viewing of adjustments.
7. Select channels and input levels to adjust, then select **up and/or down arrows**.

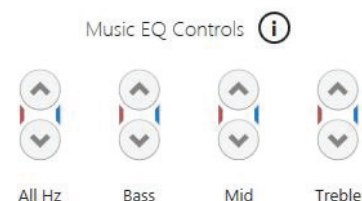

QuickFit Music Controls

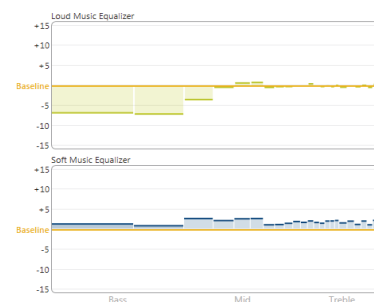

Loud and Soft Music Equalizers